

# TEST REPORT

**BioLogic Resources, LLC**  
10260 SW Nimbus Ave., Suite M11  
Portland, OR 97223  
Phone 503.670.1312  
Fax 503.670.7262

For: Yellow Emperor  
PO Box 2631  
Eugene, OR 97402

Received: 10.05.17  
Tested: 10.05.17  
Completed: 10.10.17

Lab #	Sample	APC CFU/g	Coliforms CFU/g	E. coli CFU/g	Salmonella	P. aeruginosa CFU/g	Yeast/Mold CFU/g	Listeria spp.	S. aureus CFU/g
YE145	Belly Mend	< 10	< 10	< 10	Negative	< 10	< 10 Yeast < 10 Mold	Negative	< 10
YE146	Pooper	< 10	< 10	< 10	Negative	< 10	< 10 Yeast < 10 Mold	Negative	< 10


**Method Reference:**

Bacteriological Analytical Manual (BAM), FDA, 8<sup>th</sup> edition online, 2001; Official Methods of Analysis of AOAC International, 18<sup>th</sup> edition, current through revision 4, 2011:

BAM Chapter 3: Aerobic Plate Count. Pour Plate method; AOAC Official Method 966.23  
BAM Chapter 4: Coliforms/E. coli. 3M Petrifilm; AOAC Official Method 991.14  
BAM Chapter 18: Yeast/Mold. Pour Plate method  
BAM Chapter 5: *Salmonella*. bioMerieux VIDAS UP (SPT); AOAC Official Method 2013.01  
BAM Chapter 10: *Listeria*. bioMerieux VIDAS LPT; AOAC Official Method 2013.10  
BAM Chapter 12: *S. aureus*. 3M Petrifilm; AOAC Official Method 2003.07

United States Pharmacopeia (USP), 29, 2006:  
Section <61>: *P. aeruginosa*

We warrant the above analysis was performed in good faith, using methods which are considered Standard Methods, or using methods previously agreed upon by the client. No other warranty is expressed or implied by this laboratory report. All sample results pertain only to the sample(s) tested. This test report shall not be reproduced except in full, without written approval of the laboratory.

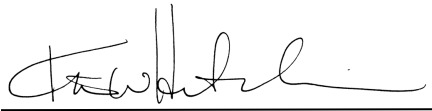

Kim W. Hutchinson, Microbiologist/Principal